




YORK AT HOME RAPID ANTIGEN TESTING PROGRAM INSTRUCTIONS

(Abbott Panbio™ package of 10 tests)

Updated: January 10, 2022

To see FAQs and distribution locations: yorku.ca/bettertogether/rapid-testing/

To upload your test result here each time: yorku.ca/go/COVID19TestResult

<p>Step 1: Gather materials</p> <ul style="list-style-type: none">• Check the expiry date on the top of this information sheet. All supplies must be discarded after that date.• You will need access to a digital camera (e.g., on a cell phone) and a smartphone, computer, or tablet with internet connection) to upload your test result.• Grab your phone or timer and a permanent marker.• Wash your hands and ensure they are completely dry before handling test supplies.• From your test package, grab 1 swab, 1 test device and 1 extraction tube and place them on a clean surface.	
<p>Step 2: Open test device</p> <ul style="list-style-type: none">• Only open the test device right before using your test, as the device is only valid for two hours after it is opened.• Make sure the testing device is laid on a flat surface with the reader facing up.	
<p>Step 3: Swab</p> <ul style="list-style-type: none">• Open the wrapper of swab like a Band-Aid®.• Only hold the plastic end of the swab and keep the cotton end facing upwards.• Insert the cotton end into your nostril (about 2 cm deep) and swirl 5 times.• Insert into second nostril and repeat.	

Step 4: Mix

- While keeping the cotton end of the swab upright, take the blue lid off the extraction tube.
- Insert the cotton end of the swab into the fluid in the extraction tube.
- Swirl the cotton tip of the swab in the fluid for 5 rotations.

**Step 5: Snap & Secure**

- Slightly raise the swab and squeeze the cotton tip to release the fluid down into the tube.
- With the swab raised, find an indentation on the plastic handle, and snap off the handle at the indentation.
- Set aside the snapped off piece to be disposed in the garbage.
- Secure blue cap on top of extraction tube with the remaining swab inside.


**Step 6: Drop**

- Twist off the bottom white lid on the extraction tube.
- Hold the tube over the round well of the test device.
- Put 5 drops into the well.
- Put the white lid back on the extraction tube.




Step 7: Wait

- Set a timer for 15 minutes.
- Wash and dry your hands.
- The test device reader window will turn completely pink – this is normal.
- Over the 15 minutes, the pink will fade to white and only the red result bar(s) will remain.

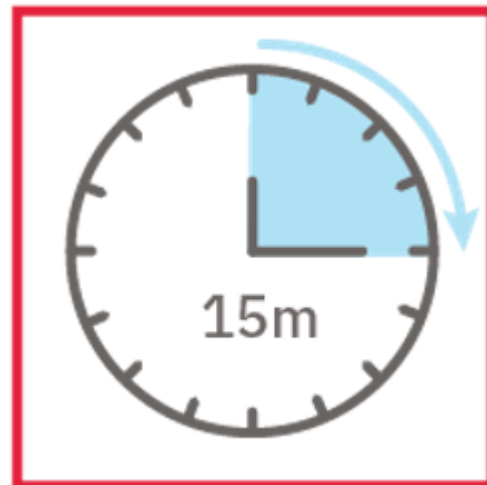


Step 8: Read your result

- After 15 minutes, check the device for your test result.
 - **Negative:** one red line at “C” level
 - **Presumed Positive:** two red lines
 - **Inconclusive:** No line or one red line at “T” level

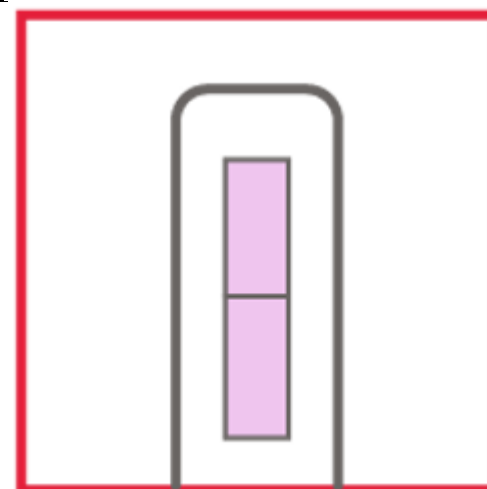

Step 9: Report your result

- With a permanent marker clearly write the test date (only date and month, e.g., Sep 25) on the test device.
- Take a photo of your test device clearly showing the result and date.
- Go to the COVID-19 Test Result form: <https://yorku.ca/go/COVID19TestResult>
- Complete the form and upload the photo of your test device.

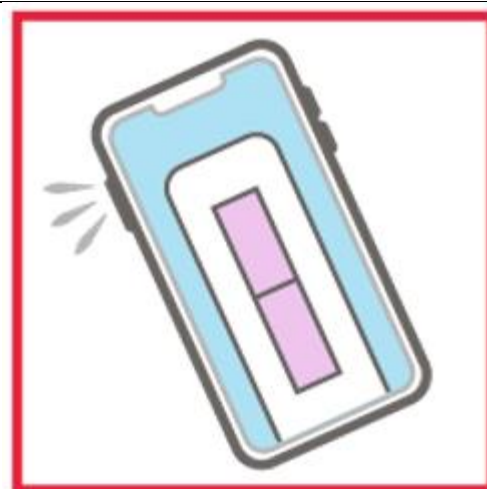

Step 10: Disposal

- Dispose of all used test materials (device, extraction tube, wrappers and swab handle). Consult with your local municipality on how to dispose of these items.

**Step 11: Next steps:****If a negative result:**

- After uploading your result to this form, no further action is required.
- If you are subject to mandatory testing, please visit the FAQs section of the [Rapid Antigen Testing Program | YU Better Together \(yorku.ca\)](#) website for more information.

If a preliminary positive result:

- Start self-isolation immediately.
- Complete the YU Screen form: <https://yorku.ubixhealth.com/> and answer **YES** on the question: “In the last 10 days, have you tested positive on a rapid antigen test or home-based self-testing kit?” A designated staff member will follow-up with you to ensure public health guidelines are being followed.

If an inconclusive result:

- Repeat the test
- If the second test result is inconclusive, please email us at: YORK-YU Testing Clinic YORK-YUTestingClinic@yuoffice.yorku.ca

Should you feel unwell at any time, complete the YU Screen to indicate that you have symptoms. If you fail the screening, stay home, do not come to campus or attend the York Testing Clinic. Self-isolate and a case and contact coordinator will be in touch with you.

York University COVID-19 testing program staff will review all uploaded results. You will be contacted should more information be required.

Resources:

- YU COVID-19 Test Result form: <https://yorku.ca/go/COVID19TestResult>
- YU rapid antigen testing program (including pick-up location information): yorku.ca/bettertogether/rapid-testing/
- YU Screen: <https://yorku.ubixhealth.com/>
- YU Better Together FAQs on mandatory testing: <https://www.yorku.ca/bettertogether/faqs/#Testing>

- Ontario COVID-19 Testing Locations: <https://covid-19.ontario.ca/assessment-centre-locations>
- Ontario Health Self-Collection training video: <https://www.youtube.com/watch?v=RXIEz1YZFDQ>
- Ontario Health Self-Collection training instructions: https://www.ontariohealth.ca/sites/ontariohealth/files/2021-03/How_to_Nasal_Swab_2-pager.pdf
- Ontario Health Considerations for Antigen Point-of-Care Testing: https://www.health.gov.on.ca/en/pro/programs/publichealth/coronavirus/docs/Antigen_Screening_Guidance_2021-03-5.pdf